


# Flexible Overtube™ Detachable Type


# sumius

For comfort repeat insertion of the endoscope during the examination and treatment.

## Structure

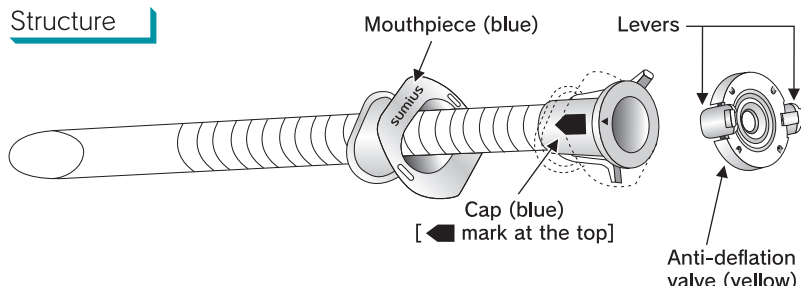

## Structure of anti-deflation valve

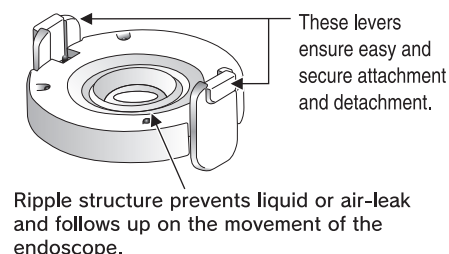

## Features

### 1. Detachable anti-deflation valve allows easy retrieval of large specimens and foreign substances.

### 2. Provides superior airtightness. (Anti-deflation valve)

- Ripple structure enables the synchronization with endoscope and maintains airtightness.

### 3. Maintains endoscopic pathway and flexibility.

- Metallic coil reinforcement prevents the bending and narrowing of the overtube while keeping its flexibility.

### 4. One touch handling.

- Mounting and demounting of the mouthpiece to the overtube.
- Fitting and detaching of the detachable valve.

Model Number	Anti-deflation valve	O.D. (mm)	I.D. (mm)	Tube length (mm)	Applicable endoscope O.D. (mm)	Packaging
MD-48718	Attached. (detachable type)	18	15	195	9~12	1 / case

Sterilized with Ethylene Oxide Gas

**Caution**

- Please do not use this product for purposes other than applicable intended use.
- This product must be used by or under the direct supervision of a doctor experienced in its use.
- Please read the attached document thoroughly before using this product.
- This product is disposable and intended for single-use. Do not re-use.
- Be sure to check the expiration date located on a sticker attached to the product.

- Please do not use the product after expiration date.
- Please do not use the product which has a damaged part in the package or a leak. Such product may be polluted.
- Store the product in a dry and cool place avoiding direct sunlight.
- If you have any question about this product, please contact our sales agency.

\*Specifications are subject to change without notice to improve technical standard, appearance or others for the product.

**SUMITOMO BAKELITE CO., LTD.**

Medical Products Business Division

2-5-8, Higashi-Shinagawa, Shinagawa-ku, Tokyo, 140-0002, JAPAN

tel: 81-3-5462-4824 fax: 81-3-5462-4894

<http://www.sumibe.co.jp>

Not FDA approved

2012.11.1000 L:AD487:12K01