

For CRYO-S ELECTRIC II

CRYOPROBES CATALOGUE


■ Gynaecology



how to identify the probe?

order/ref. no.

0 4 1024 0

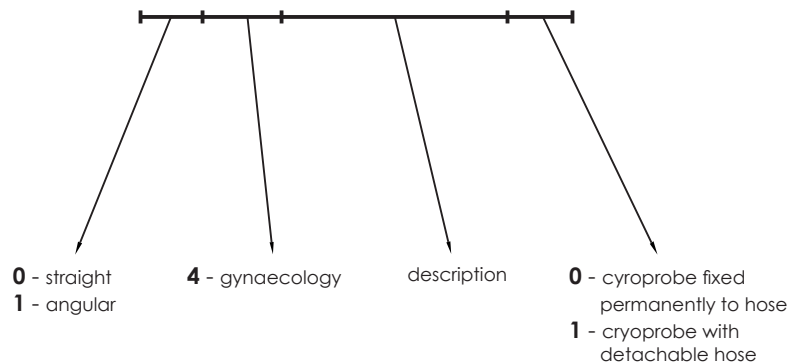

In care for highest quality of offered products and services, Metrum Cryoflex reserves the right to modify products and introduce changes in descriptions.

Gynaecology

Cryosurgery is a recognized treatment procedure using low temperatures to destroy lesions. Cryo-S Electric II unit allows controlled freezing on a measured area and depth.

Thanks to the special shapes of the cryoprobes, full scope of freezing within the area of vagina, fornix, cervix and cavum uteri is possible, as well as within the area of labia, anus and any part of the skin. For gynecology, cryosurgery is commonly used in treatment of: ectopia, polypoid and papillomatous hyperplasia of the mucous membrane, metaplasia foci with negative onco-histological result. Properly carried out freezing is painless and leaves no scars after healing. A procedure including two-time freezing lasts about 15 minutes. 4-6 weeks after the procedure there are no signs of lesions.

Cryo-S Electric II


order/ref. no. **CSEII**



State-of-the-art cryosurgical device manufactured by Metrum CryoFlex is the next generation of apparatus used by many experts in the field since 1992. Cryo-S Electric is controlled by microprocessor, and all the parameters are displayed on LCD screen.

- Mode selection, cleaning the probe and freezing can be performed automatically using footswitch keeping the site of a procedure under sterile conditions or touch screen.
- Electronic communication (chip system) between main device and connected cryoprobe. Unit recognizes optimal operating parameters and auto-configures to probe characteristics.
- Gas pressure and flow are set automatically, manual adjustment is not necessary.
- Contains automatic probe cleaning system and two freezing modes (continuous and discontinuous).
- Temperature of the probe, pressure, gas flow, time of the procedure are displayed during freezing.

CERVICAL CRYOPROBES

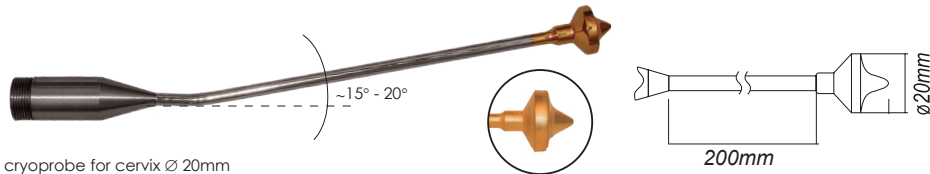

cryoprobe for cervix Ø 20mm


type: G-10/20/RF

reusable, detachable

order/ref. no. 0 4 1020 1

reusable, fixed

order/ref. no. 0 4 1020 0

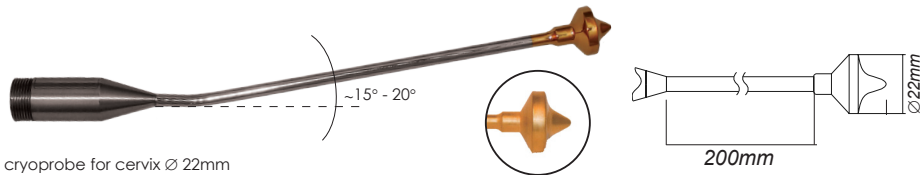

cryoprobe for cervix Ø 22mm


type: G-10/22/RF

reusable, detachable

order/ref. no. 0 4 1022 1

reusable, fixed

order/ref. no. 0 4 1022 0

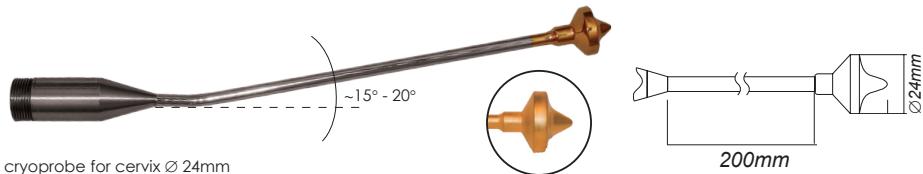

cryoprobe for cervix Ø 24mm


type: G-10/24/RF

reusable, detachable

order/ref. no. 0 4 1024 1

reusable, fixed

order/ref. no. 0 4 1024 0

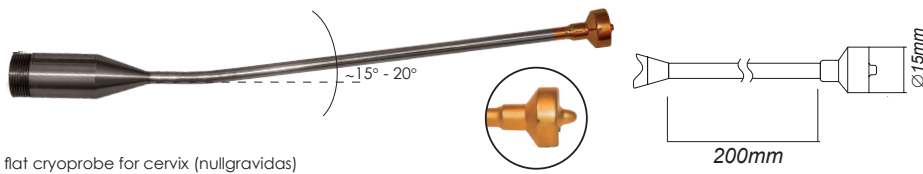

flat cryoprobe for cervix (nullgravidas)

type: G-60/RF

reusable, detachable

order/ref. no. 0 4 6015 1

reusable, fixed

order/ref. no. 0 4 6015 0

CERVICAL CRYOPROBES



multipurpose cryoprobe for additional freezing of irregular lesions

type: G-20/RF


reusable, detachable

order/ref. no. 0 4 2000 1

reusable, fixed

order/ref. no. 0 4 2000 0

ENDOMETRIUM CRYOPROBES

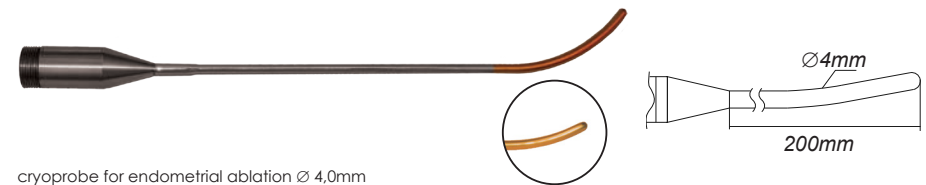

cryoprobe for endometrial ablation Ø 4,0mm


type: G-40/4/RF

reusable, detachable

order/ref. no. 0 4 4004 1

reusable, fixed

order/ref. no. 0 4 4004 0

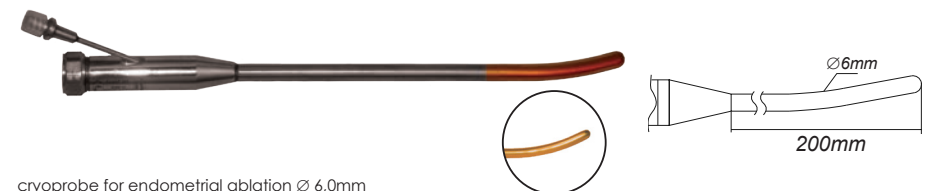

cryoprobe for endometrial ablation Ø 6,0mm

type: G-40/6/RF


reusable, detachable

order/ref. no. 0 4 4006 1

reusable, fixed

order/ref. no. 0 4 4006 0

ENDO CERVICAL CRYOPROBES

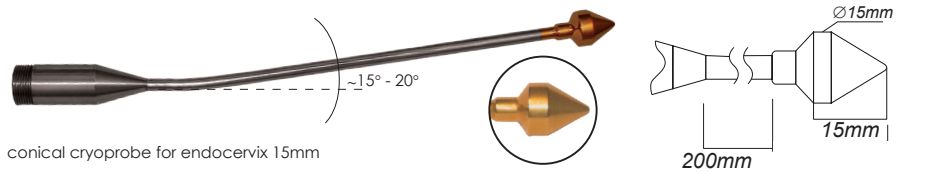


type: G-50/15/RF

reusable, detachable

order/ref. no. 0 4 5015 1

reusable, fixed

order/ref. no. 0 4 5015 0

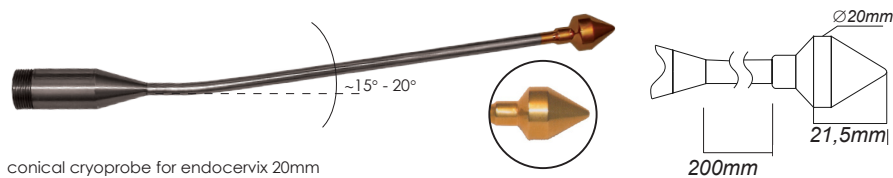


type: G-50/20/RF

reusable, detachable

order/ref. no. 0 4 5020 1

reusable, fixed

order/ref. no. 0 4 5020 0

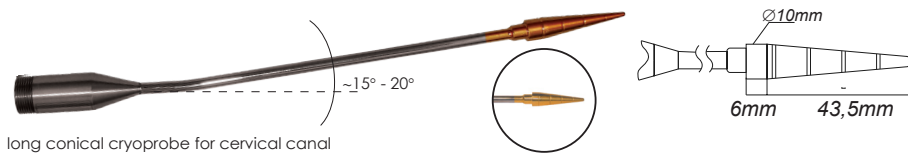


type: GH-2/RF

reusable, detachable

order/ref. no. 0 4 8002 1

reusable, fixed

order/ref. no. 0 4 8002 0

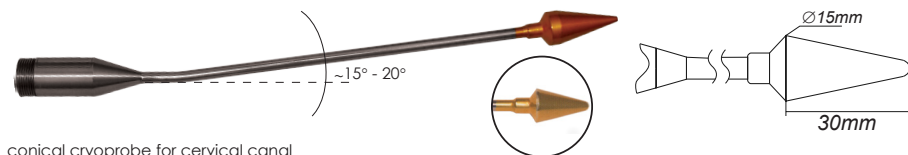


type: GH-3/RF

reusable, detachable

order/ref. no. 0 4 8003 1

reusable, fixed

order/ref. no. 0 4 8003 0

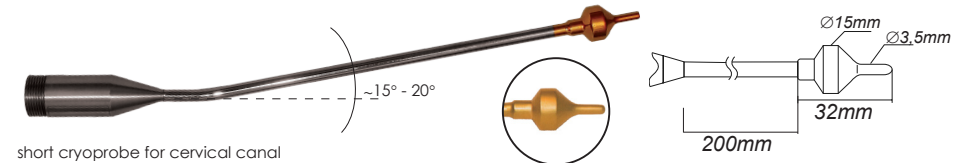ENDO CERVICAL CRYOPROBES
& CERVICAL CANAL CRYOPROBES


type: G-70/RF

reusable, detachable

order/ref. no. 0 4 7015 1

reusable, fixed

order/ref. no. 0 4 7015 0

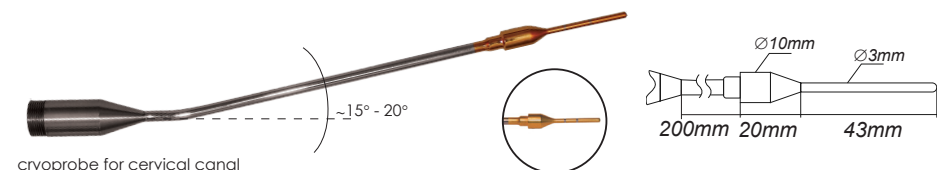


type: GH-1/RF

reusable, detachable

order/ref. no. 0 4 8001 1

reusable, fixed

order/ref. no. 0 4 8001 0

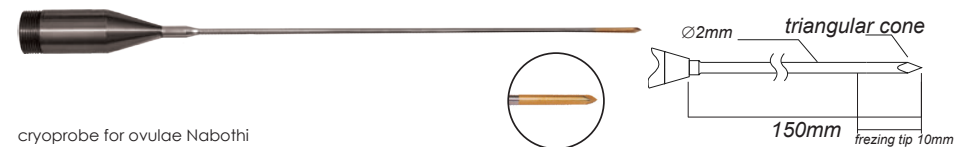


type: G-30/2/RF

reusable, detachable

order/ref. no. 0 4 3002 1

reusable, fixed

order/ref. no. 0 4 3002 0

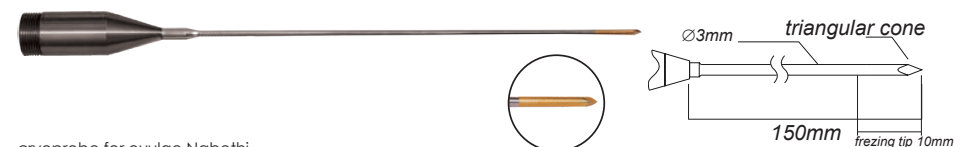

type: G-30/RF

reusable, detachable

order/ref. no. 0 4 3003 1

reusable, fixed

order/ref. no. 0 4 3003 0

SPRAY (OPEN) TYPE CRYOPROBE



multipurpose spray - type cryoprobe for external and internal application, eg: condylomata acuminata

type: **SG-10**

reusable, detachable

order/ref. no. **0 4 0000 1**



spray type cryoprobe set



multipurpose spray - type cryoprobe

type: **SD-10**

reusable, detachable

order/ref. no. **0 1 0000 1**



spray type cryoprobe set

Developed by Metrum CryoFlex spray probe cryoprobe **SG-10** is only one in the world, patented, open type cryosurgical cryoprobe with fluent, adjustable freezing power.



Reusable spray type cryoprobe can be used also with:

- single use cotton sticks (for large and regular lesions)
- reusable freezing limiters shields (for regular, flat lesions)

SG-10 spray type cryoprobe can replace all contact type cryoprobes except the needle type. Spray type cryoprobe can be used more efficiently in comparison with traditional contact type (even liquid nitrogen probes) because of:

- possibility of tissue moistening with liquid N_2O - accumulation of liquid N_2O intensifies cryodestruction (1 drop of liquid N_2O is twice as powerful as 1 drop of liquid nitrogen)
- possibility of tissue's micro cracks penetration



METRUM CRYOFLEX**Manufacture**

Zielna 29 Street
05-082 Blizne Łaszczyńskiego
POLAND, EU

Headoffice

Kolejowa 16A Street
05-092 Łomianki
POLAND, EU

Tel: +48 22 33 13 750
+48 22 33 13830
Fax: +48 22 33 13 766

export@metrum.com.pl
www.metrum.com.pl