





Designed for exceptional
kink-resistant performance
and non-ovaling effect. 

DuraSheath

Kink-resistant introducer sheath

 Rethink the basics

Supporting complex vascular access

DuraSheath

Kink-resistant introducer sheath

From easy insertion to smooth navigation to superb visibility:
DuraSheath is designed to combine optimal performance with
excellent pushability and trackability.

While maintaining the inner diameter, DuraSheath allows
for smooth device passage and has an exceptionally small
kink-radius and non-ovaling effect. The detachable hemostatic
valve provides flexibility and versatility, giving you the ability
to focus on the intervention.



Unique shaft material



newly created for safe navigation
while ensuring minimum vessel trauma.


Dilator with snap-fit

secures the dilator for safe insertion.

Detachable hemostatic valve

provides you options.

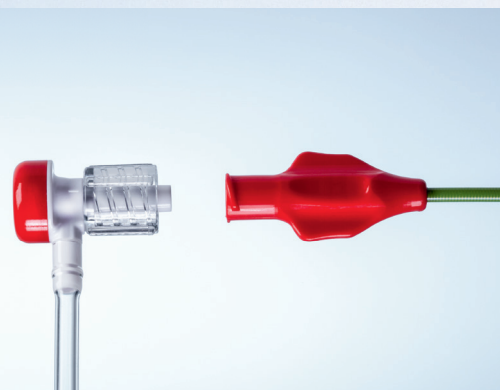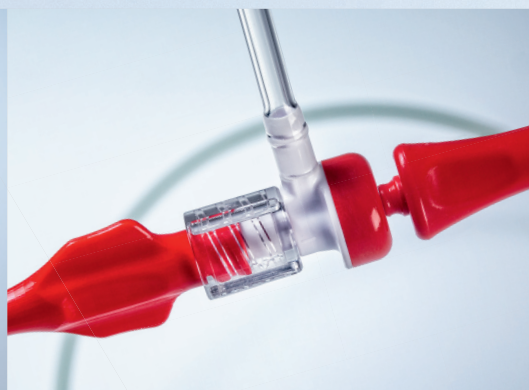

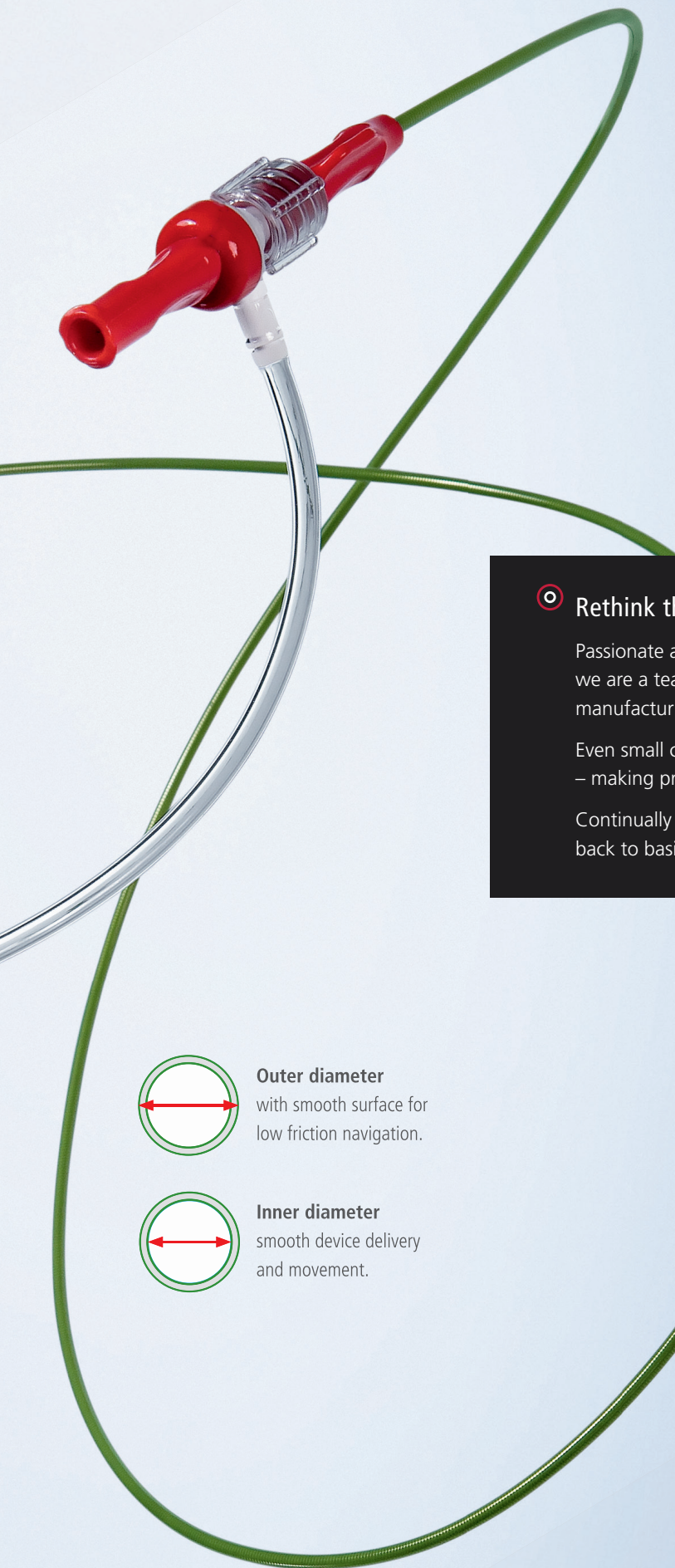


Hydrophobic coating
allowing for effortless & durable gliding
properties, no matter if wet or dry.

🎯 Rethink the basics


Passionate about bringing innovative ideas and medical device solutions to life, we are a team building on nearly 20 years of experience in the design and manufacturing of medical device solutions for partners around the world.

Even small changes in product design and quality can make a clear difference – making procedures easier and simpler with better outcomes.

Continually striving for solutions to the most basic problems, we don't go back to basics, we rethink the basics.

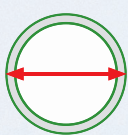


Outer diameter
with smooth surface for
low friction navigation.

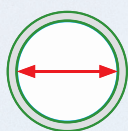


Inner diameter
smooth device delivery
and movement.



Atraumatic tip
designed for smooth and easy insertion.

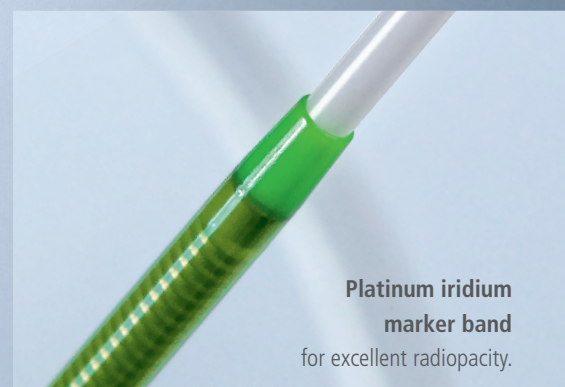

**Platinum iridium
marker band**
for excellent radiopacity.

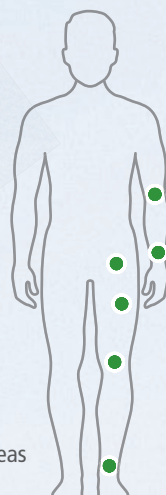

Access all areas

Technical details

Introducer Sheath

Material – outer	Stainless steel reinforced polymer
Material – inner	PTFE liner
Sheath coating	Hydrophobic
Sheath tip	Smooth taper with radiopaque platinum iridium marker band
French size	4, 5, 6, 7, 8 F
Length	15, 45, 60, 90 cm
Sidearm/hemostatic valve	Detachable, with color coded 3-way-stopcock

Ordering information

Catalogue Number	Diameter (French)	Length (cm)	Inner Diameter Minimum (mm/inch)	Dilator Guidewire Compatibility max.	Units/ Case
FG-04015	4 F	15 cm	1.60 mm / 0.063"	0.035" & 0.018"	5
FG-04045	4 F	45 cm	1.60 mm / 0.063"	0.035"	5
FG-04060	4 F	60 cm	1.60 mm / 0.063"	0.035"	5
FG-04090	4 F	90 cm	1.60 mm / 0.063"	0.035"	5
FG-05015	5 F	15 cm	1.92 mm / 0.075"	0.035" & 0.018"	5
FG-05045	5 F	45 cm	1.92 mm / 0.075"	0.035"	5
FG-05060	5 F	60 cm	1.92 mm / 0.075"	0.035"	5
FG-05090	5 F	90 cm	1.92 mm / 0.075"	0.035"	5
FG-06015	6 F	15 cm	2.25 mm / 0.088"	0.035" & 0.018"	5
FG-06045	6 F	45 cm	2.25 mm / 0.088"	0.035"	5
FG-06060	6 F	60 cm	2.25 mm / 0.088"	0.035"	5
FG-06090	6 F	90 cm	2.25 mm / 0.088"	0.035"	5
FG-07045	7 F	45 cm	2.55 mm / 0.100"	0.038"	5
FG-07060	7 F	60 cm	2.55 mm / 0.100"	0.038"	5
FG-07090	7 F	90 cm	2.55 mm / 0.100"	0.038"	5
FG-08045	8 F	45 cm	2.88 mm / 0.113"	0.038"	5
FG-08060	8 F	60 cm	2.88 mm / 0.113"	0.038"	5
FG-08090	8 F	90 cm	2.88 mm / 0.113"	0.038"	5

15 cm versions contain two dilators.
Additional sizes to become available.



A product of
Contract Medical International GmbH
Lauensteiner Str. 37
01277 Dresden • Germany
Tel. +49 351 213 88 88
www.contract-medical.com

